

FloPoint® Elite

Uroflow System

Uroflowmetry Made Easier



The FloPoint Elite Uroflow System

The FloPoint Elite Uroflow System measures and records a patient's pattern of urinary flow, allowing physicians to quickly assess the mechanics of flow and to diagnose abnormalities in voiding patterns.

The FloPoint Elite is convenient, easy to use and utilizes digital wireless technology to optimize clinical efficiency.

Portable, lightweight, and durable, the FloPoint Elite installs in seconds. Set up in a standard toilet, with no awkward chair, it frees up exam room space, helps put patients at ease and minimizes clean-up for staff.

- Helps you assess flow mechanics and diagnose voiding abnormalities (e.g., obstructed or restricted flow)
- Quickly records and charts urinary flow patterns, peak flow, voided volume, and other key measurements
- Provides readout, print-out & EMR option via HIPAA compliant ScanPoint® Technology
- Convenient online calibration
- CPT Code 51741 is approved for reimbursement by Medicare

FloPoint® Elite Uroflow System


Uroflowmetry Made Easier

FloPoint Elite Technology


The advanced technology in the FloPoint Elite Uroflow System provides an accurate and convenient method to assess bladder function. The FloSensor records and transmits patient voiding data wirelessly to the ScanPoint® Remote. This data can help determine, for example, if the patient is obstructed and requires treatment. ScanPoint technology also enables exam data storage and provides convenient exam reports.

Provides key measurements to assess bladder function

- Peak flow and average flow
- Time to Peak flow
- Voided volume
- Flow time and void time
- Flow pattern

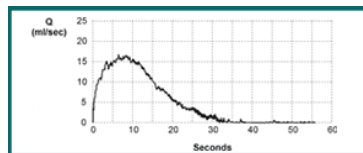

Normal Flow Pattern

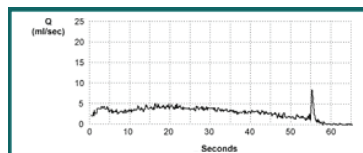

Bladder Outlet Obstruction

FloPoint Elite features:

- Quick to install, easy to clean
- Convenient for patients
- Digital flow rate
- Voice annotation
- Wireless data transfer

The FloSensor is lightweight and designed to fit in a standard toilet bowl during an exam. Between uses, it stows discreetly in the FloCharger which also recharges the FloSensor battery.

FloPoint Elite Includes:

- FloSensor and FloCharger
- ScanPoint Remote, lanyard & activation tool
- ScanPoint Docking Station
- ScanPoint In-service CD
- Funnels – 40 pack
- FloPoint Elite Manual & Quick Reference Cards
- Calibration Pack – 1 year



Sets up in a standard toilet.



Includes convenient funnels.

Dimensions:

- FloCharger – 27.5 cm high, 19.5 cm wide; weight: 1.1 kg
- FloSensor – with handles in “up” position: 26 cm high, 17 cm wide; weight: 1.05 kg
- ScanPoint Remote – 17.5 cm long, 6.25 cm wide; weight: 160 g