

GLIDESCOPE®

AVL Single Use

When time is critical, make your go-to device the GlideScope® AVL.

Ideal for demanding OR, ED, and ICU settings, the AVL provides a consistently clear view of the airway and enhanced maneuverability, enabling fast intubations.

The AVL Single Use system gives you 6 sizes to cover a wide range of patients—from preterm to morbidly obese.


Add real-time recording and an onboard tutorial, and you'll see why GlideScope AVL is the one you go to for advanced video laryngoscopy.

AVL=Advanced Airway Views

- ▶ Large digital color monitor
- ▶ Innovative blade sizing, angulation, and camera positioning
- ▶ Reveal™ anti-fog feature, with a rapid heating profile to resist lens fogging

AVL=Advanced User Features

- ▶ Real-time snapshot and video recording directly to a USB flash drive helps confirm tube placement and facilitate teaching
- ▶ Onboard *GlideScope 4-Step Technique* tutorial makes AVL easy to use, learn, and teach
- ▶ Ergonomic design and rugged construction for multiple clinical applications
- ▶ 6 sterile, single-use Stat sizes available, designed for preterm to morbidly obese patients

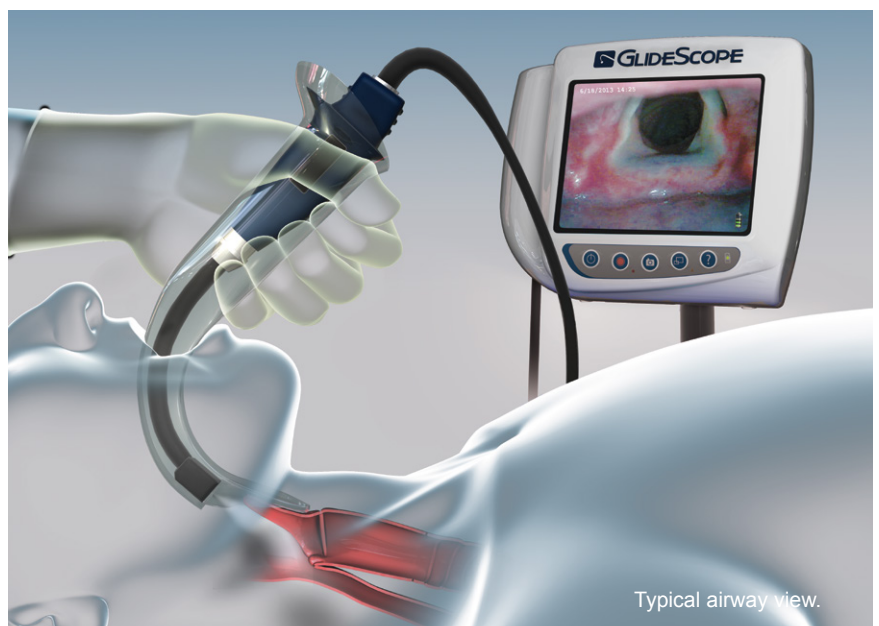

Typical airway view.











Applications

- First-use intubations, replacing direct laryngoscopy
- Normal or restricted oropharyngeal views
- Cormack-Lehane Grades I–IV laryngeal views
- Trauma airways—excellent when dealing with blood and secretions in the airway
- Nasal tracheal intubation
- Video-guided tube exchange in the ICU

System Components

- Digital, large color monitor
- Reusable video baton with digital camera
- Single-use Stat, choice of 6 sizes (box of 10)
- GlideRite® Rigid Stylet
- Protective video baton cradle
- Choice of mobile stand or IV pole mount
- Operations & Maintenance Manual and Quick Reference Card
- Instructional In-service video
- 1-year warranty

GlideScope® Products and Systems

SINGLE-USE SYSTEM				
VIDEO BATON 1-2	GVL 0	GVL 1	GVL 2	GVL 2.5
For GVL 0 - 2.5 Stats	< 1.5 kg	1.5 - 3.8 kg	1.8 - 10 kg	10 kg - 28 kg
				
VIDEO BATON 3-4	GVL 3		GVL 4	
For GVL 3 - 4 Stats	10 kg - Adult		40 kg - Morbidly Obese	
				
GLIDESCOPE PREMIUM CART			AVL PORTABLE STAND	
				

Note: System configurations vary.

800.331.2313 | verathon.com



verathon.com



Manufacturer:
Verathon Medical (Canada) ULC
 2227 Douglas Road
 Burnaby, BC V5C 5A9, Canada
 Tel: 604.439.3009
 Fax: 604.439.3039



EC Representative:
Verathon Medical (Europe) B.V.
 Willem Fenengastraat 13
 1096 BL Amsterdam, The Netherlands
 Tel: +31 (0) 20 210 30 91
 Fax: +31 (0) 20 210 30 92